



# **alpha-L-Fucosidase Microplate Assay Kit User Manual**

**Catalog # CAK1204**

(Version 1.3A)

Detection and Quantification of alpha-L-Fucosidase (AFU) Activity in  
Serum, Plasma, Tissue extracts, Cell lysate, Cell culture media and  
Other biological fluids Samples.

**For research use only. Not for diagnostic or therapeutic procedures.**

I. INTRODUCTION.....	2
II. KIT COMPONENTS.....	3
III. MATERIALS REQUIRED BUT NOT PROVIDED.....	3
IV. SAMPLE PREPARATION.....	4
V. ASSAY PROCEDURE.....	5
VI. CALCULATION.....	6
VII. TYPICAL DATA.....	7
VIII. TECHNICAL SUPPORT.....	7
IX. NOTES.....	7

## I. INTRODUCTION

alpha-L-Fucosidase (AFU) is an enzyme coded by the FUCA1 gene in humans and catalyzes the breakdown of L-Fucose. A genetic deficiency in this enzyme results in a neurovisceral storage disease, fucosidosis, which is characterized by the accumulation of fucose. Low serum activity of fucosidase has also been linked to ovarian carcinoma. Elevated fucosidase serum activity has been observed in patients with diabetes, hyperthyroidism, cirrhosis, and hepatitis. Increased activity has been associated with lung, breast, stomach, ovary, uterus, and liver carcinomas.

alpha-L-Fucosidase Microplate Assay Kit is based on the cleavage of 4-nitrophenol from the synthetic substrate. Nitrophenol becomes intensely colored after addition of the stop reagent. The increase in absorbance at 405 nm after addition of the stop reagent is directly proportional to the enzyme activity.

## II. KIT COMPONENTS

Component	Volume	Storage
96-Well Microplate	1 plate	
Assay Buffer	30 ml x 4	4 °C
Substrate	Powder x 1	-20 °C
Reaction Buffer	10 ml x 1	4 °C
Stop Solution	10 ml x 1	4 °C
Standard	Powder x 1	4 °C
Positive Control	0.1 ml x 1	-20 °C
Plate Adhesive Strips	3 Strips	
Technical Manual	1 Manual	

### Note:

**Substrate:** add 8 ml Reaction Buffer to dissolve before use.

**Standard:** add 1 ml Reaction Buffer to dissolve before use, then add 30 µl into 970 µl Reaction Buffer, mix; the concentration will be 300 µmol/L.

## III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Microplate reader to read absorbance at 405 nm
2. Distilled water
3. Pipettor, multi-channel pipettor
4. Pipette tips
5. Mortar
6. Ice
7. Centrifuge
8. Timer

#### IV. SAMPLE PREPARATION

##### 1. For cell and bacteria samples

Collect cell or bacteria into centrifuge tube, discard the supernatant after centrifugation, add 1 ml Assay buffer for  $5 \times 10^6$  cell or bacteria, sonicate (with power 20%, sonicate 3s, interval 10s, repeat 30 times); centrifuged at 8000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

##### 2. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay buffer on ice, centrifuged at 8000g 4 °C for 10 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

##### 3. For liquid samples

Detect directly.

## V. ASSAY PROCEDURE

Warm all reagents to room temperature before use.

Add following reagents into the microplate:

Reagent	Sample	Control	Standard	Blank	Positive Control
Sample	20 $\mu$ l	--	--	--	--
Distilled water	--	20 $\mu$ l	--	100 $\mu$ l	--
Positive Control	--	--	--	--	20 $\mu$ l
Standard	--	--	100 $\mu$ l	--	--
Substrate	80 $\mu$ l	80 $\mu$ l	--	--	80 $\mu$ l
Mix, put it in the oven, 37 °C for 30 minutes.					
Stop Solution	100 $\mu$ l	100 $\mu$ l	100 $\mu$ l	100 $\mu$ l	100 $\mu$ l
Mix, record absorbance measured at 405 nm.					

### Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range. If the enzyme activity is lower, please add more sample into the reaction system; or increase the reaction time; if the enzyme activity is higher, please dilute the sample, or decrease the reaction time.
- 3) Reagents must be added step by step, can not be mixed and added together.

## VI. CALCULATION

**Unit Definition:** One unit of AFU activity is defined as the enzyme generates 1  $\mu\text{mol}$  of p-nitrophenol per minute.

### 1. According to the protein concentration of sample

$$\begin{aligned}\text{AFU (U/mg)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / (C_{\text{Protein}} \\ &\quad \times V_{\text{Sample}}) / T \\ &= 0.05 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / C_{\text{Protein}}\end{aligned}$$

### 2. According to the weight of sample

$$\begin{aligned}\text{AFU (U/g)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / (V_{\text{Sample}} \times \\ &\quad W / V_{\text{Assay}}) / T \\ &= 0.05 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / W\end{aligned}$$

### 3. According to the quantity of cells or bacteria

$$\begin{aligned}\text{AFU (U/10}^4\text{)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / \\ &\quad (V_{\text{Sample}} \times N / V_{\text{Assay}}) / T \\ &= 0.05 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / N\end{aligned}$$

$C_{\text{Protein}}$ : the protein concentration, mg/ml;

$C_{\text{Standard}}$ : the concentration of Standard, 300  $\mu\text{mol/L}$  = 0.3  $\mu\text{mol/ml}$ ;

W: the weight of sample, g;

N: the quantity of cell or bacteria,  $N \times 10^4$ ;

$V_{\text{Standard}}$ : the volume of standard, 0.1 ml;


$V_{\text{Sample}}$ : the volume of sample, 0.02 ml;

$V_{\text{Assay}}$ : the volume of Assay buffer, 1 ml;


T: the reaction time, 30 minutes.

## VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

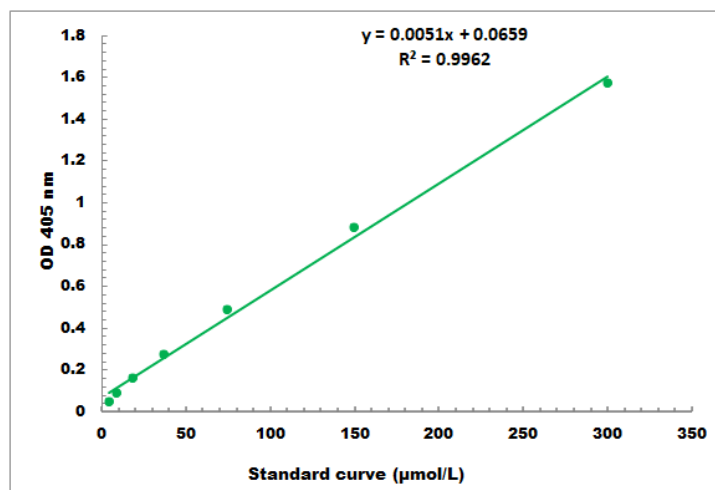

Detection Range: 3 µmol/L - 300 µmol/L

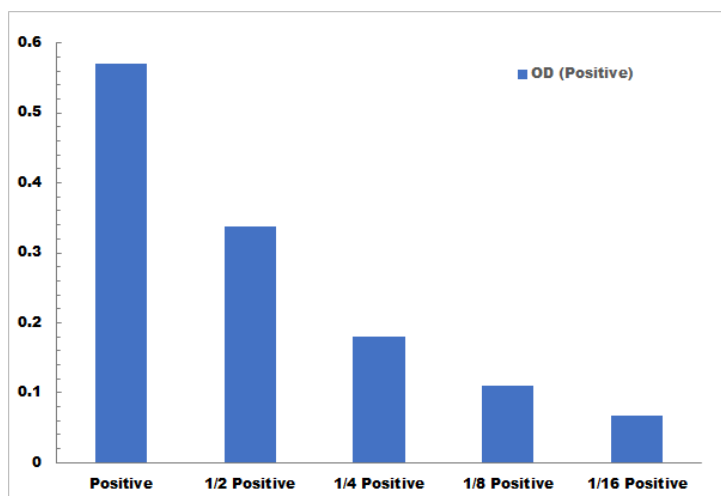

Positive Control reaction in 96-well plate assay with decreasing the concentration

## VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to [www.cohesionbio.com](http://www.cohesionbio.com) or contact us at [techsupport@cohesionbio.com](mailto:techsupport@cohesionbio.com)

## IX. NOTES