



Diamine Oxidase Activity Fluorometric Microplate Assay Kit User Manual

Catalog # CAK8008

(Version 1.1A)

Detection and Quantification of Diamine Oxidase (DAO) Activity in
Urine, Serum, Plasma, Other biological fluids, Tissue extracts, Cell
lysate, Cell culture media Samples.

For research use only. Not for diagnostic or therapeutic procedures.

| | |
|---|---|
| I. INTRODUCTION..... | 2 |
| II. KIT COMPONENTS..... | 3 |
| III. MATERIALS REQUIRED BUT NOT PROVIDED..... | 4 |
| IV. SAMPLE PREPARATION..... | 4 |
| V. ASSAY PROCEDURE..... | 5 |
| VI. CALCULATION..... | 6 |
| VII. TYPICAL DATA..... | 7 |
| VIII. TECHNICAL SUPPORT..... | 7 |
| IX. NOTES..... | 7 |

I. INTRODUCTION

Diamine oxidase (DAO) is an enzyme that your body uses to break down ingested histamine. There are a wide variety of foods that contain histamine, and it is DAO's job to break this histamine down. DAO also helps with the integrity of the gut lining, protecting us from leaky gut and the functional digestive issues that can precipitate from it.

Diamine Oxidase Activity Fluorometric Microplate Assay Kit provides a simple and sensitive method for monitoring diamine oxidase activity in various samples. In this assay, diamine is oxidized by diamine oxidase, resulting in the generation of an intermediate that reacts with the probe, which can be detected fluorometrically (Ex/Em 535/587).

II. KIT COMPONENTS

| Component | Volume | Storage |
|--------------------------|------------|----------------------|
| 96-Well Black Microplate | 1 plate | |
| Assay Buffer | 30 ml x 4 | 4 °C |
| Reaction Buffer | 20 ml x 1 | 4 °C |
| Substrate | Powder x 1 | -20 °C |
| Enzyme | Powder x 1 | -20 °C |
| Probe | Powder x 1 | -20 °C, keep in dark |
| Probe Diluent | 1 ml x 1 | 4 °C |
| Standard (100 µmol/L) | 1 ml x 1 | 4 °C |
| Positive Control | Powder x 1 | -20 °C |
| Plate Adhesive Strips | 3 Strips | |
| Technical Manual | 1 Manual | |

Note:

Substrate: add 1 ml Reaction Buffer to dissolve before use, mix. Store at 4 °C. Use within one month.

Enzyme: add 1 ml Reaction Buffer to dissolve before use. Aliquot & store at -20 °C. Use within one month.

Probe: Warm Probe Diluent to RT prior to use to melt frozen Probe Diluent; then add 1 ml Probe Diluent to dissolve. Store at -20 °C, protect from light and moisture. Use within one month.

Positive Control: add 0.1 ml Assay Buffer to dissolve before use. Store at -80 °C. Use within one month.

III. MATERIALS REQUIRED BUT NOT PROVIDED

1. Fluorescence microplate reader to read fluorescence at Ex/Em = 535/587
2. Distilled water
3. Pipettor, multi-channel pipettor
4. Pipette tips
5. Mortar
6. Centrifuge
7. Timer
8. Ice

IV. SAMPLE PREPARATION

1. For cell and bacteria samples

Collect cell or bacteria into centrifuge tube, discard the supernatant after centrifugation, add 1 ml Assay buffer for 5×10^6 cell or bacteria, sonicate (with power 20%, sonicate 3s, interval 10s, repeat 30 times); centrifuged at 10000g 4 °C for 20 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

2. For tissue samples

Weigh out 0.1 g tissue, homogenize with 1 ml Assay buffer on ice, centrifuged at 10000g 4 °C for 20 minutes, take the supernatant into a new centrifuge tube and keep it on ice for detection.

3. For liquid samples

Detect directly or dilute with Assay Buffer.

V. ASSAY PROCEDURE

Warm the solution to room temperature before use.

Add following reagents into the microplate:

| Reagent | Sample | Control | Standard | Blank | Positive Control |
|--|-------------|-------------|-------------|-------------|------------------|
| Reaction Buffer | 160 μ l | 160 μ l | 170 μ l | 170 μ l | 160 μ l |
| Sample | 10 μ l | -- | -- | -- | -- |
| Distilled water | -- | 10 μ l | -- | 10 μ l | -- |
| Standard | -- | -- | 10 μ l | -- | -- |
| Positive Control | -- | -- | -- | -- | 10 μ l |
| Substrate | 10 μ l | 10 μ l | -- | -- | 10 μ l |
| Probe | 10 μ l | 10 μ l | 10 μ l | 10 μ l | 10 μ l |
| Enzyme | 10 μ l | 10 μ l | 10 μ l | 10 μ l | 10 μ l |
| Mix, put it in the oven, 37 °C for 10 minutes, protected from light, record fluorescence measured at Ex/Em = 535/587 nm. | | | | | |

Note:

- 1) Perform 2-fold serial dilutions of the top standards to make the standard curve.
- 2) For unknown samples, we recommend doing a pilot experiment & testing several doses to ensure the readings are within the standard curve range. If the enzyme activity is lower, please add more sample into the reaction system; or increase the reaction time; if the enzyme activity is higher, please dilute the sample, or decrease the reaction time.
- 3) Reagents must be added step by step, can not be mixed and added together.

VI. CALCULATION

Unit Definition: One unit of DAO is defined as the enzyme generates 1 $\mu\text{mol H}_2\text{O}_2$ per minute at pH7.2, 37 °C.

1. According to the protein concentration of sample

$$\begin{aligned}\text{DAO (U/mg)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / \\ &\quad (C_{\text{Protein}} \times V_{\text{Sample}}) / T \\ &= 0.01 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / C_{\text{Protein}}\end{aligned}$$

2. According to the weight of sample

$$\begin{aligned}\text{DAO (U/g)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / (V_{\text{Sample}} \times \\ &\quad W / V_{\text{Assay}}) / T \\ &= 0.01 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / W\end{aligned}$$

3. According to the quantity of cells or bacteria

$$\begin{aligned}\text{DAO (U/10}^4\text{)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / \\ &\quad (V_{\text{Sample}} \times N / V_{\text{Assay}}) / T \\ &= 0.01 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / N\end{aligned}$$

4. According to the volume of sample

$$\begin{aligned}\text{DAO (U/mg)} &= (C_{\text{Standard}} \times V_{\text{Standard}}) \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}}) / V_{\text{Sample}} \\ &\quad / T \\ &= 0.01 \times (\text{OD}_{\text{Sample}} - \text{OD}_{\text{Control}}) / (\text{OD}_{\text{Standard}} - \text{OD}_{\text{Blank}})\end{aligned}$$

C_{Standard} : the concentration of standard, 100 $\mu\text{mol/L}$ = 0.1 $\mu\text{mol/ml}$;

C_{Protein} : the protein concentration, mg/ml;

W: the weight of sample, g;

N: the quantity of cell or bacteria, $N \times 10^4$;

V_{Standard} : the volume of the standard, 0.01 ml;


V_{Sample} : the volume of sample, 0.01 ml;

V_{Assay} : the volume of Assay buffer, 1 ml;


T: the reaction time, 10 minutes.

VII. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

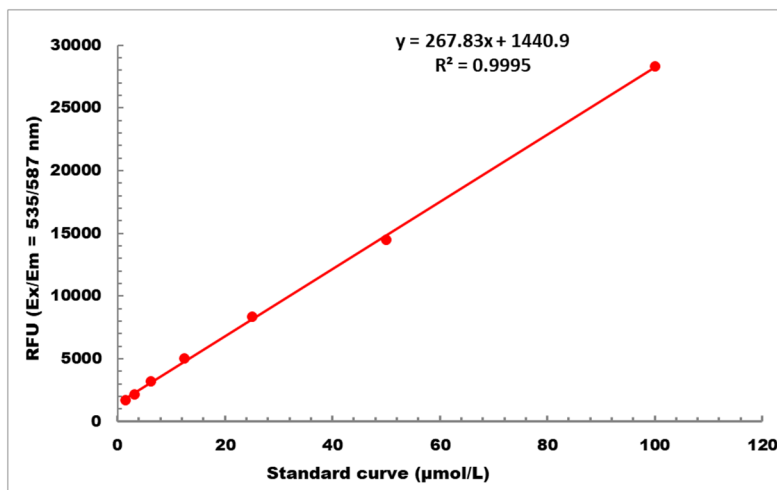

Detection Range: 1 μmol/L - 100 μmol/L

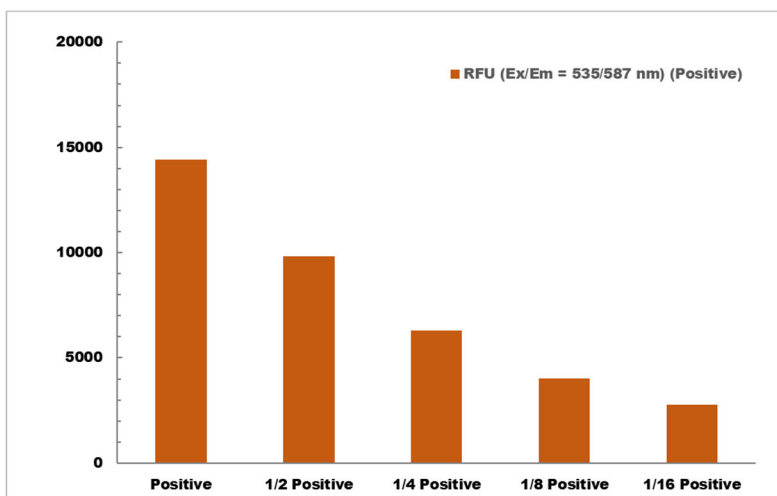

Positive Control reaction in 96-well plate assay with decreasing the concentration

VIII. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

IX. NOTES