

# Human Complement C5 ELISA Kit User Manual

Catalog # CEK1119

(Version 1.2C)

Sandwich Enzyme-Linked Immunosorbent Assay for Quantitative

Detection of Human Complement C5 Concentrations in Cell Culture

Supernatants, Serum, Plasma.

For research use only. Not for diagnostic or therapeutic procedures.



I. INTRODUCTION	2
II. ASSAY PRINCIPLES	3
III. KIT COMPONENTS	4
IV. STORAGE AND STABILITY	4
V. MATERIALS REQUIRED BUT NOT PROVIDED	5
VI. HEALTH AND SAFETY PRECAUTIONS	5
VII. REAGENT PREPARATION	6
VIII. ASSAY PROCEDURE	9
IX. ASSAY PROCEDURE SUMMARY	11
X. TYPICAL DATA	12
XI. SENSITIVITY	12
XII. SPECIFICITY	12
XIII. CROSS REACTIVITY	13
XIV. REFERENCES	13
XV. TROUBLESHOOTING GUIDE	14
XVI. TECHNICAL SUPPORT	15
XVII. NOTES	15



## I. INTRODUCTION

C5a is a protein fragment released from complement component C5. This gene is mapped to 9q33.2. The protein encoded by this gene is the fifth component of complement, which plays an important role in inflammatory and cell killing processes. This protein is comprised of alpha and beta polypeptide chains that are linked by a disulfide bridge. An activation peptide, C5a, which is an anaphylatoxin that possesses potent spasmogenic and chemotactic activity, is derived from the alpha polypeptide via cleavage with a convertase. The C5b macromolecular cleavage product can form a complex with the C6 complement component, and this complex is the basis for formation of the membrane attack complex, which includes additional complement components. Mutations in this gene cause complement component 5 deficiency, a disease where patients show a propensity for severe recurrent infections. Defects in this gene have also been linked to susceptibility to liver fibrosis and to rheumatoid arthritis.



#### II. ASSAY PRINCIPLES

The Cohesion Bioscience Human Complement C5 ELISA (Enzyme-Linked Immunosorbent Assay) kit is an in vitro enzyme-linked immunosorbent assay for the quantitative measurement of Human Complement C5 in Cell Culture Supernatants, Serum, Plasma. This assay employs an antibody specific for Human Complement C5 coated on a 96-well plate. Standards and samples are pipetted into the wells and Complement C5 present in a sample is bound to the wells by the immobilized antibody. The wells are washed and biotinylated anti-Human Complement C5 antibody is added. After washing away unbound biotinylated antibody, HRP-conjugated streptavidin is pipetted to the wells. The wells are again washed, a TMB substrate solution is added to the wells and color develops in proportion to the amount of Complement C5 bound. The Stop Solution changes the color from blue to yellow, and the intensity of the color is measured at 450 nm.



# **III. KIT COMPONENTS**

Component	Volume
96-well Plate Coated with Anti-Human Complement C5	8 wells x 12 Strips
Antibody	o wells x 12 strips
Human Complement C5 Standard	2 ng x 2
Biotin-Labeled Detection Antibody (100X)	120 μΙ
Streptavidin-HRP (100X)	120 μΙ
Standard/Sample Diluent	30 ml
Detection Antibody Diluent	12 ml
Streptavidin-HRP Diluent	12 ml
Wash Buffer (20X)	30 ml
TMB Substrate Solution	12 ml
Stop Solution	12 ml
Plate Adhesive Strips	3 Strips
Technical Manual	1 Manual

# IV. STORAGE AND STABILITY

All kit components are stable at 2 to 8 °C. Standard (recombinant protein) should be stored at -20 °C or -80 °C (recommended at -80 °C) after reconstitution. Opened Microplate Wells or reagents may be store for up to 1 month at 2 to 8 °C. Return unused wells to the pouch containing desiccant pack, reseal along entire edge. Note: the kit can be used within one year if the whole kit is stored at -20 °C. Avoid repeated freeze-thaw cycles.



# V. MATERIALS REQUIRED BUT NOT PROVIDED

- 1. Microplate reader capable of measuring absorbance at 450 nm.
- 2. Adjustable pipettes and pipette tips to deliver 2  $\mu$ l to 1 ml volumes.
- 3. Adjustable 1-25 ml pipettes for reagent preparation.
- 4. 100 ml and 1 liter graduated cylinders.
- 5. Absorbent paper.
- 6. Distilled or deionized water.
- 7. Computer and software for ELISA data analysis.
- 8. Tubes to prepare standard or sample dilutions.

#### VI. HEALTH AND SAFETY PRECAUTIONS

- 1. Reagents provided in this kit may be harmful if ingested, inhaled or absorbed through the skin. Please carefully review the MSDS for each reagent before conducting the experiment.
- 2. Stop Solution contains 2 N Sulfuric Acid ( $H_2SO_4$ ) and is an extremely corrosive agent. Please wear proper eye, hand and face protection when handling this material. When the experiment is finished, be sure to rinse the plate with copious amounts of running water to dilute the Stop Solution prior to disposing the plate.



#### VII. REAGENT PREPARATION

# 1. Sample Preparation

Store samples to be assayed within 24 hours at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.

**Cell culture supernates**: Remove particulates by centrifugation, assay immediately or aliquot and store samples at -20°C.

**Serum**: Allow the serum to clot in a serum separator tube (about 4 hours) at room temperature. Centrifuge at approximately 1000 X g for 15 minutes. Analyze the serum immediately or aliquot and store samples at -20°C.

**Plasma**: Collect plasma using heparin or EDTA as an anticoagulant. Centrifuge for 15 minutes at 1500 X g within 30 minutes of collection. Assay immediately or aliquot and store samples at -20°C.

**Cell Lysates:** Collect cells and rinse cells with PBS. Homogenize and lyse cells throughly in lysate solution. Centrifuge cell lysates at approximately 10000 X g for 5 minutes to remove debris. Aliquots of the cell lysates were removed and assayed.

**Bone Tissue:** Extract demineralized bone samples in 4 M Guanidine-HCl and protease inhibitors. Dissolve the final sample in 2 M Guanidine-HCl.

**Tissue Homogenates:** The preparation of tissue homogenates will vary depending upon tissue type. Rinse tissue with 1X PBS to remove excess blood, homogenized in 20 mL of 1X PBS and stored overnight at  $\leq$  -20 °C. After two freeze-thaw cycles were performed to break the cell membranes, the homogenates were centrifuged for 5 minutes at 5000 x g. The supernate was removed immediately and assayed. Alternatively, aliquot and store samples at  $\leq$  -20 °C.

**Note:** Some lysis buffer, such as RIPA can not be used. Some components will affect the binding.

**Urine**: Urinary samples should be cleared by centrifugation and then can be used directly without dilution. Storage at -20°C.



# 2. Human Complement C5 Standard Preparation

Dilution Buffer serves as the zero standard (0 pg/ml).

Reconstitute the lyophilized Human Complement C5 Standard by adding 0.5 ml of Standard/Sample Diluent to make the 4000 pg/ml standard stock solution. Allow solution to sit at room temperature for 5 minutes, then gently vortex to mix completely. Use within one hour of reconstituting. Two tubes of the standard (2 ng per tube) are included in each kit. Use one tube for each experiment.

Perform 2-fold serial dilutions of the top standards to make the standard curve within the range of this assay (62.5 pg/ml - 4000 pg/ml) as below. Standard/Sample

Standard Add Into 4000 pg/ml 2000 pg/ml 250 µl of the Standard (4000 pg/ml) 250 μl of the Standard/Sample Diluent 1000 pg/ml 250 µl of the Standard (2000 pg/ml) 250 μl of the Standard/Sample Diluent 500 pg/ml 250 µl of the Standard (1000 pg/ml) 250 μl of the Standard/Sample Diluent 250 μl of the Standard/Sample Diluent 250 μl of the Standard (500 pg/ml) 250 pg/ml 125 pg/ml 250 μl of the Standard (250 pg/ml) 250 μl of the Standard/Sample Diluent 62.5 pg/ml 250 µl of the Standard (125 pg/ml) 250 μl of the Standard/Sample Diluent 0 pg/ml 1 ml of the Standard/Sample Diluent

**Note:** The standard solutions are best used within 2 hours. The 4000 pg/ml standard solution should be stored at 4°C for up to 12 hours, or at -20°C for up to 48 hours. Avoid repeated freeze-thaw cycles.

# 3. Biotin-Labeled Detection Antibody Working Solution Preparation The Biotin-Labeled Detection Antibody should be diluted in 1:100 with the Detection Antibody Diluent and mixed thoroughly. The solution should be prepared no more than 2 hours prior to the experiment.

# 4. Streptavidin-HRP Working Solution Preparation



The Streptavidin-HRP should be diluted in 1:100 with the Streptavidin-HRP Diluent and mixed thoroughly. The solution should be prepared no more than 1 hour prior to the experiment.

# 5. Wash Buffer Working Solution Preparation

Pour entire contents (30 ml) of the Wash Buffer Concentrate into a clean 1,000 ml graduated cylinder. Bring final volume to 600 ml with glass-distilled or deionized water (1:20).



#### VIII. ASSAY PROCEDURE

The Streptavidin-HRP Working Solution and TMB Substrate Solution must be kept warm at 37°C for 30 minutes before use. When diluting samples and reagents, they must be mixed completely and evenly. Standard detection curve should be prepared for each experiment. The user will decide sample dilution fold by crude estimation of protein amount in samples.

- 1. Add 100 µl of each standard and sample into appropriate wells.
- 2. Cover well and incubate for 90 minutes at room temperature or over night at 4°C with gentle shaking.
- 3. Remove the cover, discard the solution and wash plate 3 times with Wash Buffer Working Solution, and each time let Wash Buffer Working Solution stay in the wells for 1 2 minutes. Blot the plate onto paper towels or other absorbent material. Do NOT let the wells completely dry at any time.
- 4. Add 100  $\mu$ l of Biotin-Labeled Detection Antibody Working Solution into each well and incubate the plate at 37°C for 60 minutes.
- 5. Wash plate 3 times with Wash Buffer Working Solution, and each time let Wash Buffer Working Solution stay in the wells for 1 2 minutes. Discard the Wash Buffer Working Solution and blot the plate onto paper towels or other absorbent material.
- 6. Add 100  $\mu$ l of Streptavidin-HRP Working Solution into each well and incubate the plate at 37°C for 45 minutes.
- 7. Wash plate 5 times with Wash Buffer Working Solution, and each time let wash buffer stay in the wells for 1 2 minutes. Discard the wash buffer and blot the plate onto paper towels or other absorbent material.
- 8. Add 100  $\mu$ l of TMB Substrate Solution into each well and incubate plate at 37°C in dark for 10-20 minutes.
- 9. Add 100  $\mu$ l of Stop Solution into each well. The color changes into yellow immediately.



10. Read the O.D. absorbance at 450nm in a microplate reader within 30 minutes after adding the Stop Solution.

For calculation, (the relative O.D.450) = (the O.D.450 of each well) - (the O.D.450 of Zero well). The standard curve can be plotted as the relative O.D.450 of each standard solution (Y) vs. the respective concentration of the standard solution (X). The concentration of the samples can be interpolated from the standard curve.

**Note:** If the samples measured were diluted, multiply the dilution factor to the concentrations from interpolation to obtain the concentration before dilution.




# IX. ASSAY PROCEDURE SUMMARY

M	Prepare all reagents, samples and standards
M	• Add 100 μl Standard or Sample
M	Wash plate 3 times with Wash Buffer Working Solution
M	• Add 100 μl Biotin-Labeled Detection Antibody Working Solution
M	Wash plate 3 times with Wash Buffer Working Solution
M	• Add 100 μl Streptavidin-HRP Working Solution
M	Wash plate 5 times with Wash Buffer Working Solution
M	• Add 100 μl TMB Substrate Solution
M	• Add 100 μl Stop Solution
	• Read the plate at 450nm



# X. TYPICAL DATA

The standard curve is for demonstration only. A standard curve must be run with each assay.

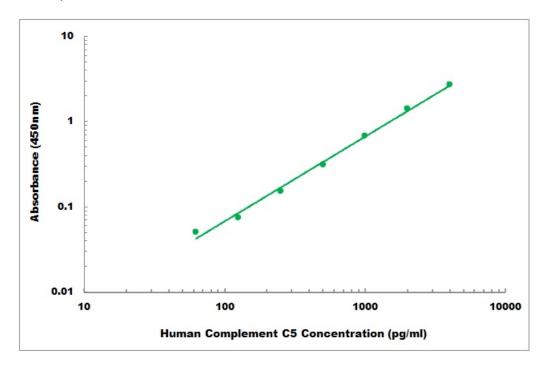

# XI. SENSITIVITY

The minimum detectable dose of Human Complement C5 is typically less than 30 pg/ml.

# XII. SPECIFICITY

The Human Complement C5 ELISA Kit allows for the detection and quantification of endogenous levels of natural and/or recombinant Human Complement C5 proteins within the range of 62.5 pg/ml - 4000 pg/ml.



# XIII. CROSS REACTIVITY

No detectable cross-reactivity with other relevant proteins.

## XIV. REFERENCES

- 1. Delgado-Cervino, E., Fontan, G., Lopez-Trascara, M. C5 complement deficiency in a Spanish family: molecular characterization of the double mutation responsible for the defect. Molec. Immun. 42: 105-111, 2005.
- 2. Pfarr, N., Prawitt, D., Kirschfink, M., Schroff, C., Knuf, M., Habermehl, P., Mannhardt, W., Zepp, F., Fairbrother, W. G., Loos, M., Burge, C. B., Pohlenz, J. Linking C5 deficiency to an exonic splicer enhancer mutation. J. Immun. 174: 4172-4177, 2005. Note: Erratum: J. Immun. 182: 5152 only, 2009.



# XV. TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
High signal and background in	Insufficient washing	Increase number of washes
all wells		Increase time of soaking
		between in wash
	Too much Streptavidin-HRP	Check dilution, titration
	Incubation time too long	Reduce incubation time
	Development time too long	Decrease the incubation
		time before the stop solution
		is added
No signal	Reagent added in incorrect	Review protocol
	order, or incorrectly prepared	
	Standard has gone bad (If	Check the condition of
	there is a signal in the sample	stored standard
	wells)	
	Assay was conducted from an	Reagents allows to come to
	incorrect starting point	20 - 30 °C before performing
		assay
Too much signal-whole plate	Insufficient washing-unbound	Increase number of washes
turned uniformly blue	Streptavidin-HRP remaining	Carefully
	Too much Streptavidin-HRP	Check dilution
	Plate sealer or reservoir	Use fresh plate sealer and
	reused, resulting in presence of	reagent reservoir for each
	residual Streptavidin-HRP	step
Standard curve achieved but	Plate not developed long	Increase substrate solution
poor discrimination between	enough	incubation time
point	Improper calculation of	Check dilution, make new
	standard curve dilution	standard curve
No signal when a signal is	Sample matrix is masking	More diluted sample
expected, but standard curve	detection	Recommended
looks fine		
Samples are reading too high,	Samples contain protein levels	Dilute samples and run
but standard curve is fine	above assay range	Again
Edge effect	Uneven temperature around	Avoid incubating plate in
	work surface	areas where environmental
		conditions vary
		Use plate sealer



# XVI. TECHNICAL SUPPORT

For troubleshooting, information or assistance, please go online to www.cohesionbio.com or contact us at techsupport@cohesionbio.com

# **COHESION BIOSCIENCES LIMITED**

FLAT32 ADVENTURES COURT
12 NEWPORT AVENUE
LONDON, E14 2DN, UK

Website: www.cohesionbio.com

Email: order@cohesionbio.com

techsupport@cohesionbio.com

custom@cohesionbio.com

XVII. NOTES